


THERMOTRAK II

Control up to 32 Thermotron environmental test chambers from a single computer using ThermoTrak II software. This free software provides a private test chamber network and is your central hub for all chamber programming, operating, data acquisition, and display functions.

ThermoTrak II software can log data as fast as every two seconds and is only limited by the amount of hard drive space on your computer. Controlling environmental test chambers is more efficient than ever before, thanks to ThermoTrak II's innovative capabilities that simplify lab management and improve performance by saving time, increasing productivity, and cutting costs.

This software was designed to help improve the accuracy, speed, and quality of your test results. It displays all the real-time variables your chamber controllers provide, offering full-color graphs of the setpoints, process variables, and monitor channels. ThermoTrak II's automated features monitor and report various chamber events and alarms, so users are aware of what is happening while multiple tests are taking place. The software was designed to help improve the accuracy, speed, and quality of your test results.

ThermoTrak II has a Windows®-based interface, which makes it easy to install and use, and is compatible with current and most legacy Thermotron controllers. ThermoTrak II is free and can be downloaded from the Thermotron website. You can set up your own private network of environmental chambers by following the steps below.

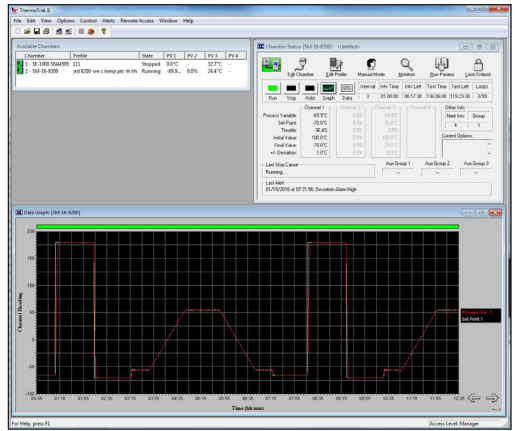

Controller Requirements	
Controller Type	Computer Interface
2800 / 4800	RS-232 RS-485 GPIB TCP/IP
3200	
3800	
6850	
7800	
7825	
8200	
8800/8825	
Oven	

1. Open your web browser and navigate to: thermotron.com/service-2/controller-software-updates-downloads/
2. Download the software by right-clicking on the file name and choose 'Save Target As...'
3. Open the file using zip-capable software such as WinZip.

For further instructions, view the "readme" file using Microsoft Notepad.

If you have questions or would like more information, please contact a Thermotron technical advisor or your local sales representative.

For more than 55 years, Thermotron has provided quality environmental test equipment. We've worked to establish a trusted reputation among our peers, and when people hear the name *Thermotron*, they have confidence in the testing of their own product. We've been building our name since 1962; now it's your turn.

**QUALITY. TRUST.
CONFIDENCE.**

— BUILD YOURS WITH A —

THERMOTRON.

THERMOTRON.COM

US: 291 Kollen Park Drive, Holland, Michigan 49423 | P: (616) 393-4580 | F: (616) 392-5643 | info@thermotron.com

UK: Winch Rd., Kent Science Park, Sittingbourne, Kent, ME9 8EF England | P: 01795 436333 | F: 01795 436777 | sales@thermotron.co.uk